


AMT NutraGlide®

NG/NJ Nasal Feeding Tubes

**NEW
SIZES**

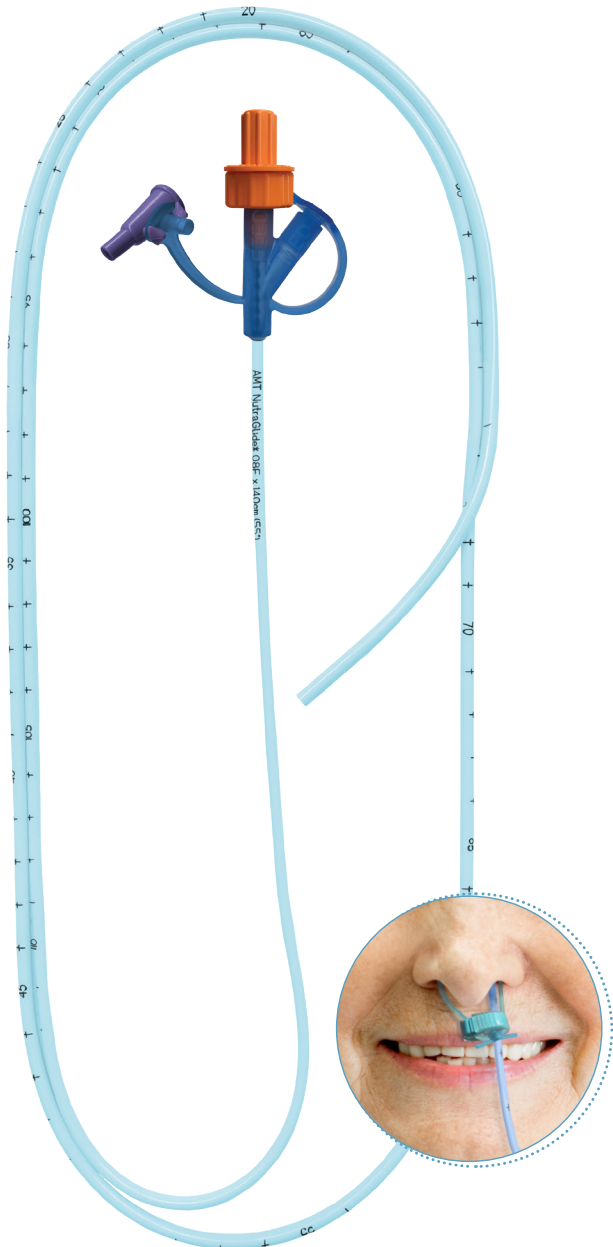

Designed to Reduce Clogging

A unique increasing internal diameter and open distal tip is designed to reduce clogging.

Decrease Replacement Tubes

When compared to the NutraGlide®, a straight tube is more than twice as difficult to expel a clog formation.¹

Improve Comfort

Smaller proximal diameter of the nasal tube within the nasal passage “may improve patient comfort and likely decrease the incidence of sinusitis.”¹

Safe + Secure

When used together, the Bridle Pro® with NutraGlide® improve nutrition, promote recovery, and decrease costs. Bridle Every Tube, Every Time!

Ready Glide® Technology Eliminates Prewetting for Stylet Removal!



Innovating. Educating. Changing Lives.™


800 869 7382 | CS@AppliedMedical.net | ICS@AppliedMedical.net

AMT NutraGlide®

Changing the Standard of Care for Nasal Feeding

NutraGlide® is Designed to Reduce Clogging

The inner and outer diameters of the NutraGlide® Nasal Feeding Tube gradually increase from the proximal entry port to the distal tip, which allows for improved flow and reduces the risk of clog formation. The reverse taper tube design also makes it easier to remove a clog, if one should form.¹ The 5F NutraGlide® transitions from 5F at the proximal entry port to 6F at the distal tip; similarly, the 6F transitions from 6F to 8F, the 8F transitions to 10F, and the 10F NutraGlide® transitions to 12F.

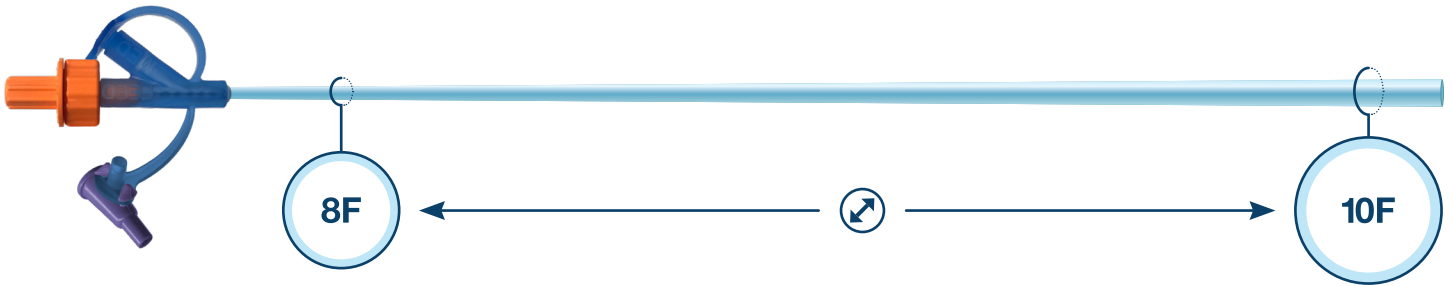

Ready Glide® Technology

All sizes feature AMT's proprietary Ready Glide® coating, which eliminates the need for water activation, and is intended to provide an easier placement and removal of the stylet.

NutraGlide® | Order Information

FRENCH SIZE	LENGTH	LEGACY		ENFIT®	
		STYLET	NO STYLET	STYLET	NO STYLET
5F	38 cm (15")	NFT-05038	NFT-05038-N	NFT-05038-I	NFT-05038-N-I
	55 cm (21")	NFT-05055	NFT-05055-N	NFT-05055-I	NFT-05055-N-I
	90 cm (35")	NFT-05090	NFT-05090-N	NFT-05090-I	NFT-05090-N-I
6F	55 cm (21")	NFT-06055	NFT-06055-N	NFT-06055-I	NFT-06055-N-I
	90 cm (35")	NFT-06090	NFT-06090-N	NFT-06090-I	NFT-06090-N-I
8F	90 cm (35")	NFT-08090	NFT-08090-N	NFT-08090-I	NFT-08090-N-I
	110 cm (43")	NFT-08110	NFT-08110-N	NFT-08110-I	NFT-08110-N-I
	140 cm (55")	NFT-08140	NFT-08140-N	NFT-08140-I	NFT-08140-N-I
10F	110 cm (43")	NFT-10110	NFT-10110-N	NFT-10110-I	NFT-10110-N-I
	140 cm (55")	NFT-10140	NFT-10140-N	NFT-10140-I	NFT-10140-N-I

Add the prefix "S" to order sterile Nasal Feeding Tubes, e.g. SNFT-05038. | NutraGlide® is sold 10/box. The NutraGlide® is guidewire compatible throughout the entire length of the device.

For a complete list of products, visit us at www.AppliedMedical.net
800 869 7382 | CS@AppliedMedical.net | ICS@AppliedMedical.net



www.AppliedMedical.net



Applied Medical Technology, Inc.
 8006 Katherine Boulevard
 Brecksville, OH 44141
www.AppliedMedical.net
 An ISO 13485 Registered Company
 MR Safe (with stylet removed)
 Made In USA



© 2025 Applied Medical Technology, Inc. **C-4739-G**

Patents: <http://www.AppliedMedical.net/Patents/> | Trademarks: <http://www.AppliedMedical.net/Trademarks/>

ENFit® is a registered trademark of GEDSA. | 1. Kirn, D. S., Whitman, W. L., Hisel, R. D., & Price, J. (2022, June). Analysis of a Novel Nasoenteral Tube Design. In The American Society of Mechanical Engineers. <https://asmedigitalcollection.asme.org/medicaldevices/article-abstract/doi/10.1115/1.4053876/1136724/Analysis-of-a-Novel-Nasoenteral-Tube-Design?redirectedFrom=fulltext>