



MiniACE® Button Cecostomy + Appendicostomy

MiniACE® Button Low Profile Antegrade Continence Enema Button

MiniACE® Enema Button | Order Information

Length (cm)	ORDER NUMBERS					
	10F		12F		14F	
	LEGACY	ENFIT®	LEGACY	ENFIT®	LEGACY	ENFIT®
0.8	ACE-1008	ACE-1008-I	ACE-1208	ACE-1208-I	ACE-1408	ACE-1408-I
1.0	ACE-1010	ACE-1010-I	ACE-1210	ACE-1210-I	ACE-1410	ACE-1410-I
1.2	ACE-1012	ACE-1012-I	ACE-1212	ACE-1212-I	ACE-1412	ACE-1412-I
1.5	ACE-1015	ACE-1015-I	ACE-1215	ACE-1215-I	ACE-1415	ACE-1415-I
1.7	ACE-1017	ACE-1017-I	ACE-1217	ACE-1217-I	ACE-1417	ACE-1417-I
2.0	ACE-1020	ACE-1020-I	ACE-1220	ACE-1220-I	ACE-1420	ACE-1420-I
2.3	ACE-1023	ACE-1023-I	ACE-1223	ACE-1223-I	ACE-1423	ACE-1423-I
2.5	ACE-1025	ACE-1025-I	ACE-1225	ACE-1225-I	ACE-1425	ACE-1425-I
2.7	ACE-1027	ACE-1027-I	ACE-1227	ACE-1227-I	ACE-1427	ACE-1427-I
3.0	ACE-1030	ACE-1030-I	ACE-1230	ACE-1230-I	ACE-1430	ACE-1430-I
3.5	ACE-1035	ACE-1035-I	ACE-1235	ACE-1235-I	ACE-1435	ACE-1435-I
4.0	ACE-1040	ACE-1040-I	ACE-1240	ACE-1240-I	ACE-1440	ACE-1440-I
4.4	ACE-1044	ACE-1044-I	ACE-1244	ACE-1244-I	ACE-1444	ACE-1444-I
5.0	ACE-1050	ACE-1050-I	ACE-1250	ACE-1250-I	ACE-1450	ACE-1450-I
5.5	ACE-1055	ACE-1055-I	ACE-1255	ACE-1255-I	ACE-1455	ACE-1455-I
6.0	ACE-1060	ACE-1060-I	ACE-1260	ACE-1260-I	ACE-1460	ACE-1460-I
6.5	ACE-1065	ACE-1065-I	ACE-1265	ACE-1265-I	ACE-1465	ACE-1465-I
7.0	ACE-1070	ACE-1070-I	ACE-1270	ACE-1270-I	ACE-1470	ACE-1470-I
8.0	ACE-1080	ACE-1080-I	ACE-1280	ACE-1280-I	ACE-1480	ACE-1480-I
9.0	ACE-1090	ACE-1090-I	ACE-1290	ACE-1290-I	ACE-1490	ACE-1490-I
10.0	ACE-10100	ACE-10100-I	ACE-12100	ACE-12100-I	ACE-14100	ACE-14100-I

MiniACE® Button comes packaged as Qty/Box 1 (sterile). | Kit contents include: device (1), irrigation set (1), 5cc syringe (1), introducer/stiffener (1), gauze, and water soluble lubricant (1). MR Conditional.

Irrigation Set Included in Kit



Complete Portfolio – Irrigation Sets

ADAPTER OPTIONS	RIGHT ANGLE CONNECTOR OPTIONS	
	12"	24"
Bolus	2-1221	2-2421
	2-1221-ISOSAF	2-2421-ISOSAF

Note: Irrigation sets are packaged as 10/box (non-sterile)

Balloon Fill Volumes

FRENCH SIZE	MIN	REC.	MAX
10F	1 ml	1.5 ml	2 ml
12F	2 ml	2.5 ml	3 ml
14F	3 ml	4 ml	5 ml

Recommended water fill volume is clearly printed on the balloon fill valve of each MiniACE® Button. Consult your healthcare provider before adjusting fill volumes.


Balloon Measuring Device | Order Number 1904

Guidewire Compatible (up to 0.038"), 10F Tubing




AMT CINCH® Tube Securement

Order Numbers: CINCH410S, CINCH417M, and CINCH930L

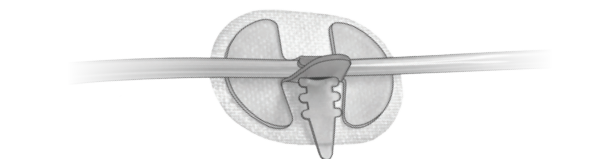

AMT Clamp™ | Order Number 4-3000

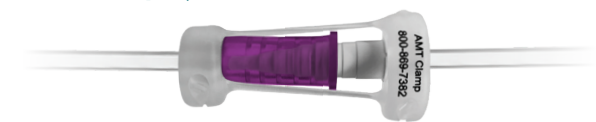

"With the MiniACE®, she empties her bowels daily, she no longer has a distended stomach, and she is pain-free!"

Audrey's Mom

AppliedMedical.net/Stories/Audreys-MiniACE-Enema-Button/



Cecostomy + Appendicostomy

Cleared for Initial Placement

Superior Leakage Control

Improved Stoma Retention

For a complete list of products, visit us at www.AppliedMedical.net
800 869 7382 | CS@AppliedMedical.net | ICS@AppliedMedical.net



MiniACE® Button



Applied Medical Technology, Inc.
8006 Katherine Boulevard
Brecksville, OH 44141
www.AppliedMedical.net
An ISO 13485 Registered Company
Cleared For Initial Placement
MR Conditional
Made in USA



© 2026 Applied Medical Technology, Inc. C-4642-F

Patents: <http://www.AppliedMedical.net/Patents/> | Trademarks: <http://www.AppliedMedical.net/Trademarks/>
Chait Trapdoor™ is a trademark of Cook® Medical | ENFit® is a registered trademark of GEDSA. | 1.Claim based on internal testing of 10F-14F devices.



Innovating. Educating. Changing Lives.™

Secure Attachment



Size Variety




Comfortable




MiniACE® Button

- Cleared for Initial Placement
- Minimizes stoma dilation associated with conversion from Chait Trapdoor™ to MiniACE®
- Cecostomy + Appendicostomy

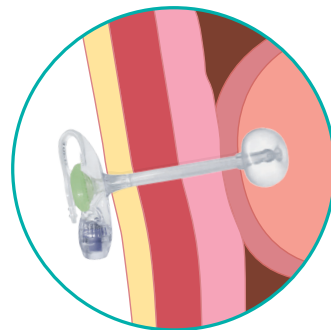

Cecostomy

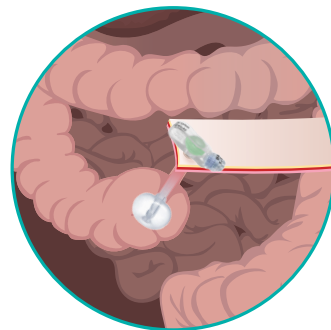


Appendicostomy

Patient Comfort First

- Medical grade silicone, gentle to the skin
- Flexible + soft, easy to clean around the button
- Exclusive AMT Balloon, designed to reduce pullouts and leakage
- Balloon bolster shields the internal tip of the MiniACE®, protecting the cecum


Reduce Leaks, Superior Retention

- AMT's MiniACE® Button is 450% more secure within the stoma when compared to the Chait Trapdoor™.¹
- Exclusive AMT Balloon, along with tapered tubing at the external bolster, reduce the risk of leaks.
- Anti-leakage valve helps prevent backflow.

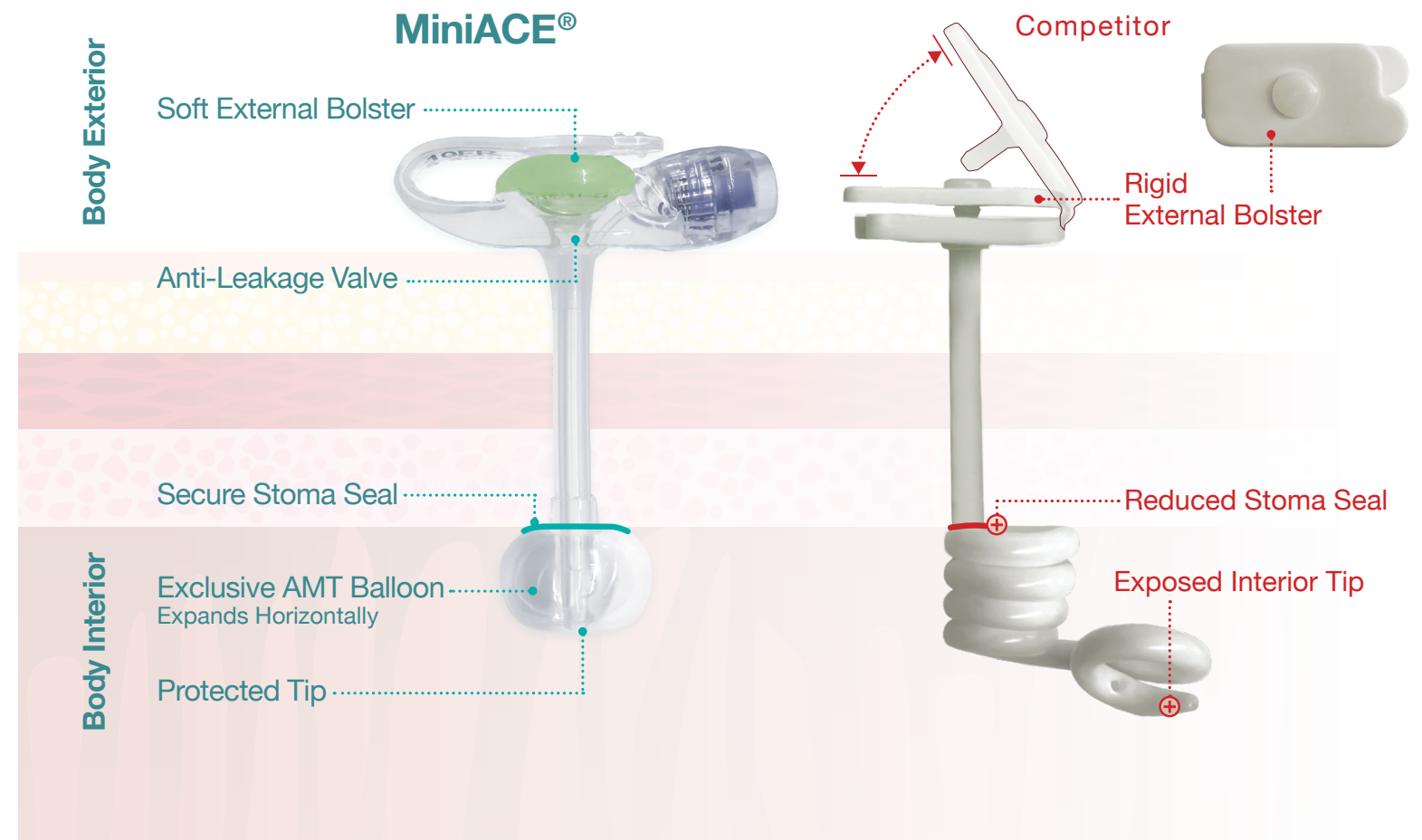


Safe and Secure Irrigation Set Attachment


- The connection between the MiniACE® Button and its irrigation set is 10x stronger than the competition.¹

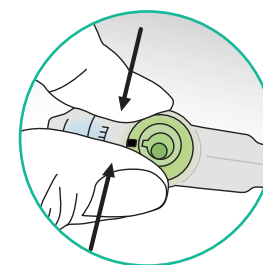


STABILIZE

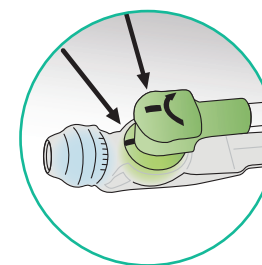

ALIGN & INSERT

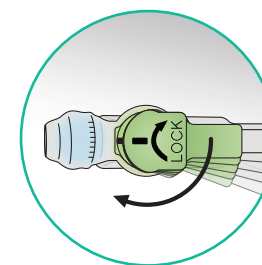

ROTATE

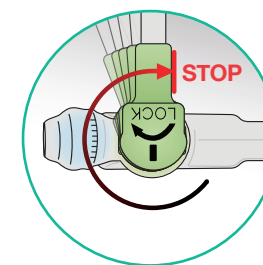

LOCK



IRRIGATE

